



BARAN FARMIN TADBIR



OLIDENT
INNOVATION

Composite REVOLUTION

RIGHT SOLUTION FOR EVERY SITUATION



New

- revolutionary OliREVO
- expert OliCo XP
- classic VIP esthetic
- complimentary:
VIP Posterior
- standard OliFlow
- **NEW OliFlow SE**
www.BFT-co.com




Made in Germany



We are proud to present the top three composite materials among their competitors. Each of them has different features, together they have one thing in common - all meet the strictest German quality standards.

German Quality Guaranteed!

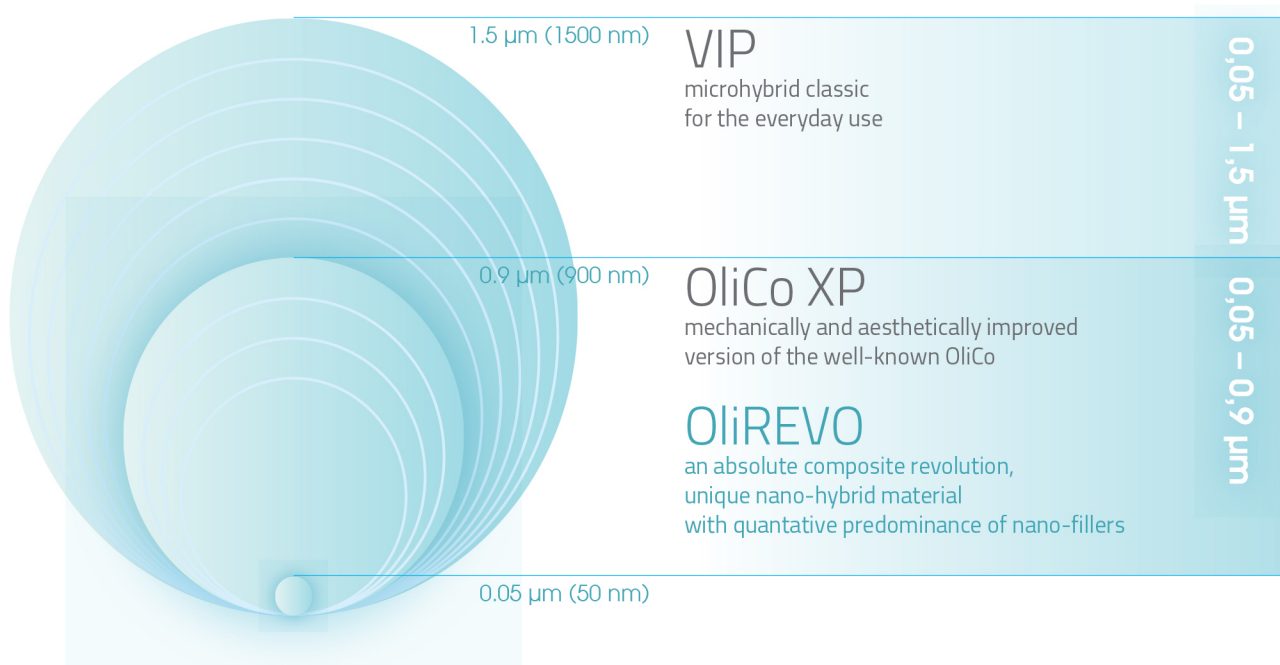

INDICATIONS	COMPOSITE						
	OliREVO	OliCo XP	VIP esthetic	VIP Posterior	OliFlow	OliFlow SE	
Class I	●	●	●	●	●	●	
Class II	●	●	●	●	●	●	
Class III	●	●	●				
Class IV	●	●	●				
Class V	●	●	●		●	●	
Minimally invasive techniques					●	●	
Primary teeth fillings	●	●	●	○	○	○	
Direct composite veneers	●	●	○				
Inlays, onlays and indirect composite veneers	●	●	○				
Baseliner					●	●	
Sealing pits and fissures					●	●	

● main indication

○ additional indication




OliREVO

First on the market nano-hybrid composite, that is so convenient and easy to form. Its consistency is called „Wax-like“ enhancing the creation of highly aesthetic restorations with mechanical properties that guarantee the long-term reliable results.


OliREVolution! Absolute composite

LONG-TERM AESTHETIC FILLINGS

Thanks to the exceptional resistance to compression **OliREVO** tolerates high loads, both on the occlusal and incisal surfaces. In addition, **OliREVO** has an unusually low solubility in water - more than seven times lower than the ISO standard! This and the low wear of **OliREVO** provides for long-term esthetic restorations.

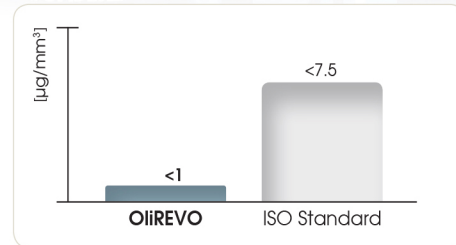

OliREVO water solubility index is more than 7 times lower than the ISO standard*
*Manufacturer internal measurements, Germany.

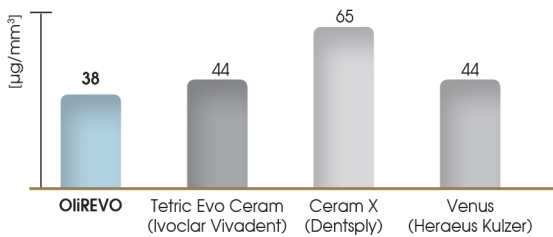

Abrasive wear of OliREVO fillings and competitive materials (200 000 cycles of abrasion ~ approximately one year)*
*Manufacturer internal measurements, Germany.

ABRASIVE WEAR RESISTANCE

The studies on the fillings surface abrasive wear show that **OliREVO** achieves one of the best results among its competitors. After 200 000 cycles the fillings surface presented 38 μm wear. This proves **OliREVO** restorations are extremely durable and also safe for the opposing teeth.

„WAX-like“ CONSISTENCY FOR DENTIST COMFORT

The German specialists who created **OliREVO** took great care to make a user friendly composite that will surely impress all users. The unique WAX-like consistency is synonymous with smooth and easy handling of **OliREVO**. The non-stick characteristics allow for efficient molding and carving even in difficult restorative surfaces.

WAX like

OliREVO - clinical picture by courtesy of Dr. Dominik Morawski, *Denticulus, Nowy Sacz/Limanowa*



Patient with diastema - initial situation.



Diastema closure with OliREVO, tooth 12,11 and 21,22..




New
New

OliCo XP

OliCo XP is a technologically improved generation of OliCo, a composite already known to the market. Along with excellent esthetic advantages, there are also impeccable mechanical properties that existing and new users will appreciate. **OliCo XP** is an ideal material for both anterior and posterior teeth.

OliCo - an eXPert
in nano-ceramic technology

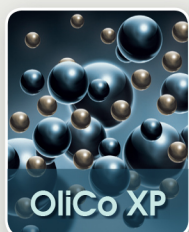


The new **OliCo XP** has more nano-particles with preserved overall filler content than OliCo.



MUCH MORE „NANO” FILLER

The new **OliCo XP** has far more nano-particles with preserved overall filler content. Modification of the structure helped to obtain much better strength parameter and lower polymerization shrinkage. **OliCo XP** fillings are more durable and abrasion resistant.

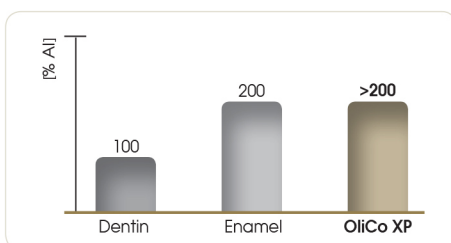

Exceptional visibility in radiographs.

FULL RANGE OF BENEFITS AND APPLICATIONS

- ✓ abrasion resistance
- ✓ low shrinkage
- ✓ pronounced „chameleon effect” for excellent aesthetics
- ✓ ease of polishing and long-lasting shine
- ✓ non-stick effect aids in ease of molding and shaping
- ✓ wide range of applications

Clinical pictures by courtesy of Dr. Barbara and Tomasz Smigiel, *Smigiel Dentistry, Katowice*



The secondary carious lesions, tooth 11.



Removed carious lesions.



The finished composite restoration.



OliArtist Brush

Forming composite with special brushes is a commonly used technique. Here we present OliArtist Brush made from the highest quality Kolinsky natural bristles and available in several sizes. They are an absolute must when a particular aesthetic shape and color of the filling is desired.

Sizes: #1 - OL 838, #2 - OL 839

Art of aesthetics



VIP esthetic


For many years **VIP esthetic** has been a favorite microparticle hybrid composite in the market. Excellent colour stability and mechanical durability was gained thanks to very low rate of water sorption.

VIP esthetic microhybrid technology classic



VIP ESTHETIC - EXCEPTIONALLY LOW WATER SORPTION


Recent research of the **VIP esthetic** properties shows a very low rate of water sorption. General ISO standards for composites require index lower than 40 $\mu\text{g}/\text{mm}^3$, **VIP esthetic** was rated 4 times lower! For this reason, VIP restorations show excellent color stability and mechanical durability.

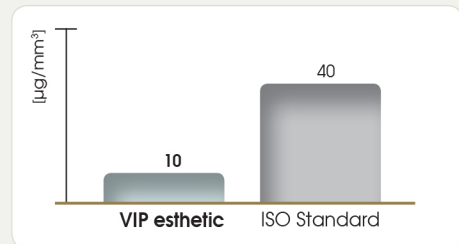

Water sorption of the **VIP esthetic** composite is 4 times lower than ISO standard*

*Manufacturer internal measurements, Germany.

VIP esthetic application - clinical pictures by courtesy of Dr. Dominik Morawski, *Denticulus, Nowy Sacz/Limanowa*

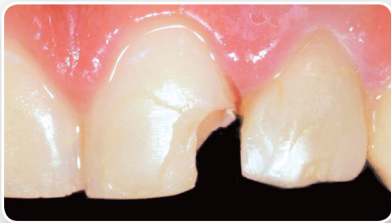

Removed carious lesions, tooth 21.



The finished composite restoration.

VIP Posterior


Due to its excellent mechanical strength and resistance to abrasion, **VIP Posterior** is intended for posterior restorations that are particularly exposed to the occlusion forces. Increased filler content (weight: 84%) makes the polymerization shrinkage significantly lower and thus stress affecting the cavity walls is minimized.

OPTIMUM FILLINGS HARDNESS

In comparison with VIP esthetic, which is a universal composite, **VIP Posterior** was specially designed for posterior fillings. The filler content was increased in order to gain higher hardness of the material but without losing its excellent consistency. **VIP Posterior** is available in 3 colors: A2, A3 and A3.5.



New
New

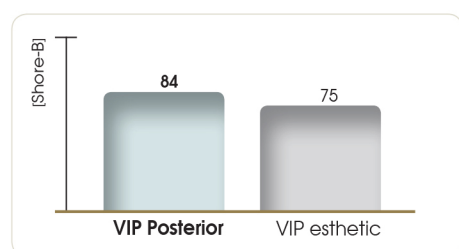

VIP Posterior hardness on the Barcol scale is much higher than the hardness of VIP esthetic*
*Manufacturer internal measurements, Germany.

Composite
for special tasks

New




OliFlow SE

Self adhesive light-curing flowable composite for wide range of indications in the restorative therapy.

OliFlow SE has self adhesives properties, so requires no etching, priming and bonding of enamel and dentine. The procedure time is much shorter, minimizes any potential mistakes and ensures postoperative sensitivity reduction.

etching
priming
bonding
filling

All in One

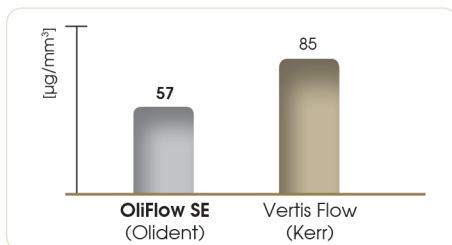

Abrasive wear of **OliFlow SE** and competitive material (200 000 cycles of abrasion - approximately 1 year)*

*Manufacturer internal measurements, Germany.

ABRASIVE WEAR RESISTANCE

Low wear of a filling material is essential for clinical success. Abrasive wear of **OliFlow SE** after 200 000 cycles presents 57µm. It guarantees durable restorations over long term use.

APPLICATION INDICATIONS FOR OliFlow SE AND OliFlow

- ✓ base lining in deep cavities
- ✓ pit and fissure sealings
- ✓ small restorations of class I and II
- ✓ all types of restorations in primary teeth
- ✓ enamel defects corrections
- ✓ blocking out and filling of undercuts
- ✓ temporary restoration corrections



fluoride
release



OliFlow

A proven and appreciated flowable microhybrid composite designed for small cavities, fissure sealing and a liner material in deep restorations. Thanks to a great consistency it adapts perfectly to the cavity walls preventing flow outside the preparation walls.

MEETS THE TWO MOST IMPORTANT PROTECTION FEATURES:

- ✓ protects the pulp against mechanical and thermal factors
- ✓ provides insulation against the irritating effects of substances contained in the fillings and the external environment



OliEtch BOND




Light cured, one component, single layer, VII generation bonding system. Created to simplify bonding procedures while maintaining excellent quality and strong bond strength. **One Layer - Zero Mistakes!**


ONE OF THE STRONGEST - LOOK AND COMPARE

Among the most popular self-etching systems, **OliEtch BOND** is distinguished by excellent bonding strength values. Only one layer of **OliEtch BOND** is needed to get strong and long-lasting bonding of the light-curing composites to enamel and dentin.*

The bond strength to enamel (MPa)*

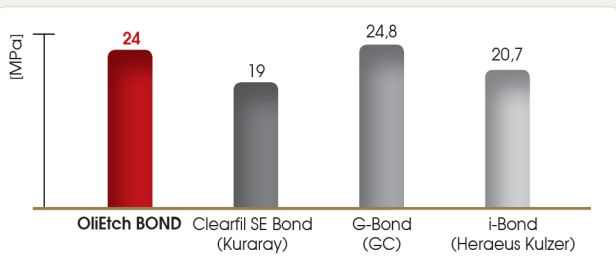

The bond strength to dentin (MPa)*

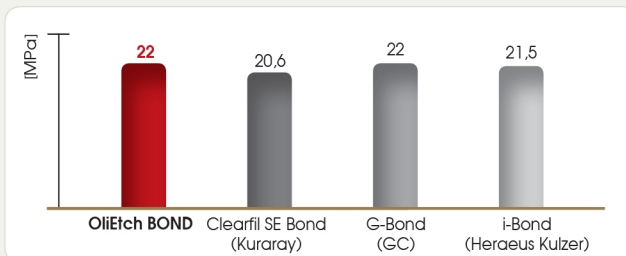


*Manufacturer internal measurements, Germany.

ONE procedure in 50 seconds!


Earlier generations of bonding systems required the use of several components / bottles and several separate procedures. With **OliEtch BOND** all the procedures: etching, conditioning and bonding take place in one stage and takes only 50 seconds!




Apply a small amount of **OliEtch BOND** on the tray.

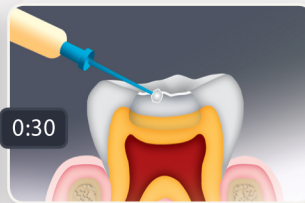

Use the applicator to thoroughly rub **OliEtch BOND** in moist enamel and dentin for 30 seconds

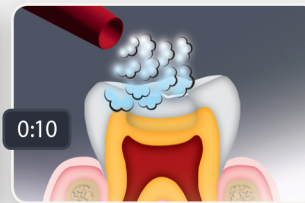

Blow air for 10 seconds to remove any excess **OliEtch BOND**.

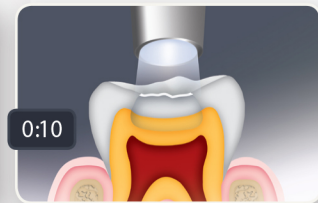

Polymerize **OliEtch BOND** with a curing- light for 10 seconds.



OliBOND / OliActivator

OliBOND is a two step Prime and Bond system that is light cured providing a strong bond of all composite materials to dentin and enamel. Combine **OliBOND** with **OliActivator** and you now have a self curing chemical bond as well as a dual cure system. Ideal for root canal and post placement as well as other dual cure materials.

Reliable in every clinical situation!



OLIDENT



BARAN FARMIN TADBIR
www.BFT-co.com

Composite REVolution

RIGHT SOLUTION FOR EVERY SITUATION

New
New

OliREVO single

1 syringe 4g in any shade:

SHADE	A1	A2	A3	A3.5	OA2	OA3	BL	B1	B2	Incisal
CODE:	OL 311	OL 312	OL 313	OL 314	OL 315	OL 322	OL 394	OL 316	OL 317	OL 321
ECO:	OL 311E	OL 312E	OL 313E	OL 314E	OL 315E	OL 322E	OL 394E	OL 316E	OL 317E	OL 321E

OliREVO set 7***

7 syringes 4g in any shades:

CODE:	OL 373
ECO:	OL 373E

***Your own composition of shades!

Only with us!

Selected shades are also available in single doses 0,3g.

New
New

OliCo XP single

1 syringe 5g in any shade:

SHADE	A1	A2	A3	A3.5	OA2	OA3	BL	B1	B2	C2	D2	D3	Incisal
CODE:	OL 011XP	OL 012XP	OL 013XP	OL 014XP	OL 015XP	OL 022XP	OL 094XP	OL 016XP	OL 017XP	OL 018XP	OL 019XP	OL 020XP	OL 021XP
ECO:	OL 011XPE	OL 012XPE	OL 013XPE	OL 014XPE	OL 015XPE	OL 022XPE	OL 094XPE	OL 016XPE	OL 017XPE	OL 018XPE	OL 019XPE	OL 020XPE	OL 021XPE

OliCo XP set 12***

12 syringes 5g in any shades:

CODE:	OL 071XP
ECO:	OL 071XPE

OliCo XP set 8***

8 syringes 5g in any shades:

CODE:	OL 076XP
ECO:	OL 076XPE

VIP esthetic single

1 syringe 5g in any shade:

SHADE	A1	A2	A3	A3.5	OA2	B2
CODE:	OL 211	OL 212	OL 213	OL 214	OL 215	OL 217
ECO:	OL 211E	OL 212E	OL 213E	OL 214E	OL 215E	OL 217E

VIP esthetic set 12***

12 syringes 5g in any shades:

CODE:	OL 271
ECO:	OL 271E

VIP esthetic set 8***

8 syringes 5g in any shades:

CODE:	OL 276
ECO:	OL 276E

OliFlow SE

1 syringe 2ml (3g)

New

SHADE	A2
CODE:	OL 058

VIP Posterior single

1 syringe 5g in any shade:

New
New

SHADE	A2	A3	A3.5
CODE:	OL 252	OL 253	OL 254
ECO:	OL 252E	OL 253E	OL 254E

OliFlow

1 syringe 2ml (3g) in any shade:

SHADE	A1	A2	A3
CODE:	OL 043	OL 042	OL 040

OliEtch BOND OliBOND

1 bottle 5ml

1 bottle 5ml

OliActivator

1 bottle 5ml

CODE:	OL 003
ECO:	OL 003E

CODE:	OL 001
ECO:	OL 001E

CODE:	OL 002
ECO:	OL 002E



Attention! Join us and „be eco“!

Take care for the environment with us!

When you buy selected products from our catalog, you can opt for paper-free packaging. In return you will receive „eco-discount“. Make your environment cleaner and save money!

Your business Partner



Join us on Facebook

tel. +98 21 77 955 687

info@BFT-co.com

www.BFT-co.com



Olident is not responsible for any printing errors.



OLIDENT

EN ISO
13485:2012
+AC:2012



www.BFT-co.com